


### TECHNICAL SPECIFICATIONS:

Characteristic	SculptIris
Lipolysis Method	Heat-assisted
Laser Type	Diode
Wavelength	1060nm
Energy	072-1.87w/cm <sup>2</sup>
Pulse Length	1s-20s
Work time	25 minutes
Applicator size	55*105mm
Treatment clamp indicating light	Red (630nm)/Green(570nm)/50mW×4
Working Temperature	5°C~30°C
Package Size	120*62*60mm & 60*46*19mm
Net Weight	75&11kg
Rated Input Power	AC110V/60Hz AC220V/50Hz
Maximum Input Power	2500W

*A simple and fast way  
to reduce fat!*

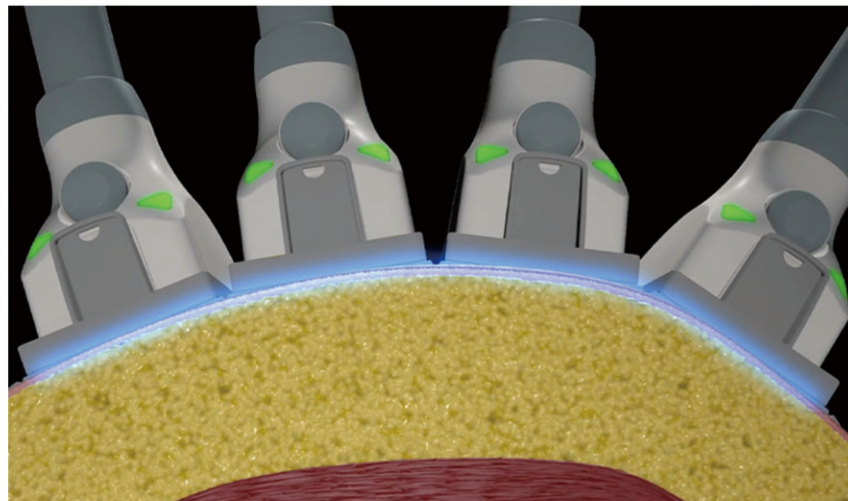

1060nm

The 1060nm wavelength has a special affinity for adipose tissue, it can effectively treat the stubborn cellulite.

# SculptIris

Non-invasive Body Slimming with  
1060nm Laser





Lumzeus®

Medical CE FDA ISO 13485

# SculpTiris


## Theory:

SculpTiris is the world's first FDA-cleared Technology laser treatment for non-invasive lipolysis of the abdomen, flanks, back, inner thighs, outer thighs and is now available for the submental (under the chin) area. It's designed to safely and effectively eliminate unwanted fat cells. The result is a permanent reduction in stubborn fat without surgery or downtime.

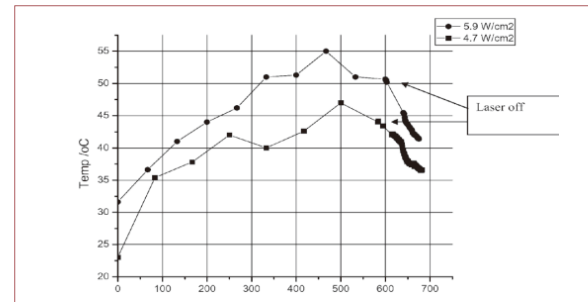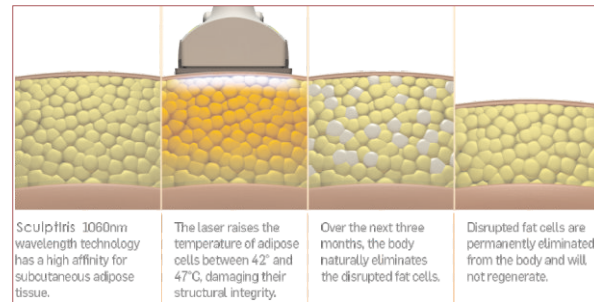

## Characteristics

- Minimal absorption in the dermis leaves the skin's surface unharmed
- Advanced Contact Cooling enhances patient comfort
- Feathering of heat spread provides natural-looking results
- Fast, 25-minute treatment per area
- Flat, non-suction applicator design for consistent and natural-looking results
- Simple and efficient attachment system
- Non-invasive customized treatments for your clients
- No surgery, no downtime


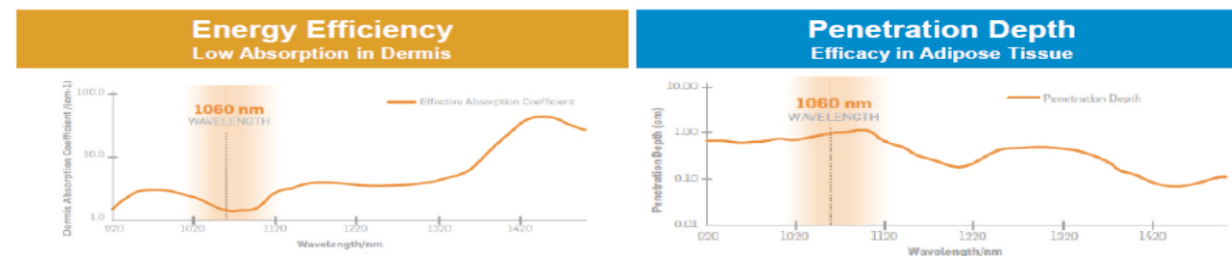


## A Smarter Way to SculpTiris

- Flat, non-suction applicator design for consistent, natural looking results
- Non-invasive customize treatment for your patients
- Simple and efficient attachment system

## The areas that can be treated:

- Under the chin
- Belly
- Love handles
- Back
- Inner and outer thighs

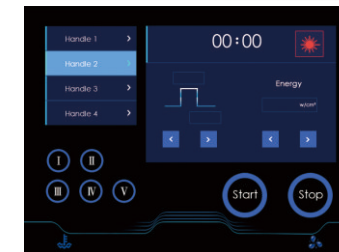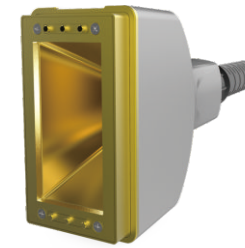

The 1060 nm wavelength's specific affinity for adipose tissue, coupled with minimal absorption in the dermis, allows SculpTiris to efficiently treat areas of troublesome fat in just 25 minutes. The body naturally eliminates the disrupted fat cells with results seen as quickly as 6 weeks and optimal results usually seen in as few as 12 weeks.

## Before and After:

